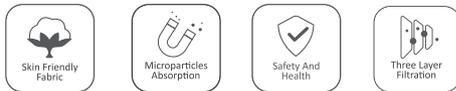
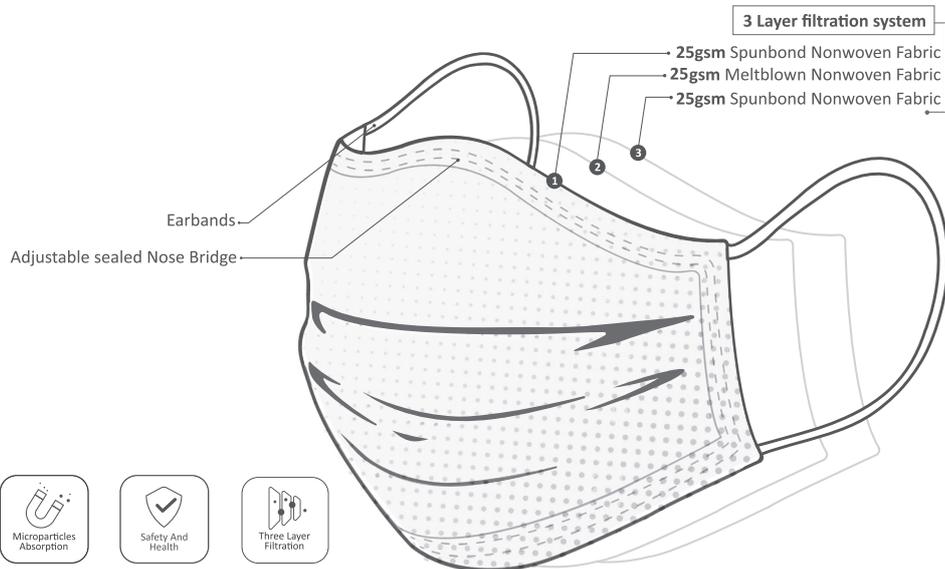


A 3ply face mask is a loose-fitting, disposable device that creates a physical barrier between the mouth and nose of a 3ply face mask is a loose-fitting, disposable device that creates a physical barrier between the mouth and nose of the wearer and potential contaminants in the immediate environment.


MANUFACTURER:

InSpe
Savanoriu ave. 178F, Vilnius, Lithuania

PRODUCT

MEDICAL FACE MASK

PRODUCT NAME:

BALTIC MASKS

CLASSIFICATION:

Type IIR

MODEL:

BM-220 (Disposable 3-ply Medical Non Sterile Face Mask)

QUALITY ASSURANCE AND CONFORMITY:

Directive 93/42/EEC

Regulation: **(EU) 2017/745**
of the European Parliament and of the Council of 5 April 2017 on
Medical Devices

EN 14683:2019+AC:2019 Annex B

Bacterial Filtration Efficiency (BFE)

EN 14683:2019+AC:2019 Annex C

Breathability (Differential Pressure)

ISO 22609:2004

Resistance against penetration by synthetic bloo

EN ISO 11737-1: 2018

Microbial cleanliness (Bioburden)

Test report: **AR-20-YL-004978-02 (2020 08 26) EUROFIN**

MATERIALS:

Outer Material	1st layer - 25gsm Spunbond Non-Wowen
Filter Layer	2nd layer - 25gsm Melt-Blown fabric
Inner Material	3rd layer - 25gsm Spunbond Non-Wowen
Nose-piece	Plastic single wire thread 3mm
Earloops	3mm Polyster and Spandex

SIZE:

Mask body **175x95mm (+/- 5mm)**

SHELF LIFE:

3 years after production date

STORAGE:

Storage under cool, clean, and dry conditions. Avoid excessive heat (-20°C to +40°C and at less than 80% relative humidity)

 Usage instruction

1. Before opening and touching the mask, wash your hands thoroughly with soap and water or alcohol/based hand sanitizer (>60% alcohol content).
2. Open the mask with both hands to the two sides.
3. Hold the mask by the ear loops. Place a loop around each ear and be sure the nose clip is upwarded
4. Press the nose clip to fit the shape of the nose
5. Pull the bottom of the mask over your mouth and chin

 Precautions

1. This product is not suitable for protection against harmful gases, steam and oily aerosols.
2. The product does not generate oxygen and is not suitable for use in anoxic environments, underwater operations, during evacuation or firefighting environments.
3. It is not recommended to use this protective mask if you do not feel a good adhesion between the face and the edge of the mask during use.
4. If found during use that the mask is damaged or the breathing resistance is too large, replace with a new one.
5. This product is for one-time use, masks should not be reused.

